Update

Reference Guide Directions for Use

PROPAQ *encore*[®]
Vital Signs Monitor

English Update for software versions 2.0X-2.1X.
Keep this update in the back pocket of your manual.

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 **PROTOCOL**[®]
SYSTEMS, INC.

Copyright © 2000, Protocol Systems, Inc.

Protocol[®], Propaq[®], Propaq Encore[®], Acuity[®], Protocol Cordless[®], and Flexible Monitoring[®] are registered trademarks and Smartcuf[™] is a trademark of Protocol Systems, Inc. Protocol Systems, Inc. is protected under various patents and patents pending.

Copyright Notice: Software in this Product is Copyright© 2000 by Protocol Systems, Inc., or its vendors. All rights are reserved. The software is protected by United States of America copyright laws and international treaty provisions applicable all over the world. Under such laws, you are licensed to use the copy of the software incorporated with this instrument as intended in the operation of the product in which it is embedded, but the software may not be copied, decompiled, reverse-engineered, disassembled or otherwise reduced to human-perceivable form. This is not a sale of the software or any copy of the software; all right, title and ownership of the software remains with Protocol or its vendors. Protocol Systems, Inc. will make available specifications necessary for interoperability of this software on request; however, users should be aware that use of Protocol hardware and software with devices or software not sold by Protocol or its authorized dealers and affiliates may lead to erroneous results and consequent danger in patient care, and may also void Protocol's warranty.

DISCLAIMERS:

Protocol Systems, Inc. cautions the reader of this manual:

- This manual may be wholly or partially subject to change without notice.
- All rights are reserved. No one is permitted to reproduce or duplicate, in any form, the whole or part of this manual without permission from Protocol Systems, Inc.
- Protocol Systems, Inc. will not be responsible for any injury to the user or other person(s) that may result from accidents during operation of the Propaq Encore.
- Protocol assumes no responsibility for usage not in accordance with this manual that results in illegal or improper use of the Propaq Encore.

Send correspondence regarding this manual and other Protocol products to:

Protocol Systems, Inc.

Customer Service
8500 SW Creekside Place
Beaverton, Oregon 97008-7107 USA

Within USA, toll free:

Phone: (800) 289-2500
Phone Technical Services: (800) 289-2501

Worldwide:

Phone: (503) 526-8500
Fax: (503) 526-4200
Fax Technical Services: (503) 526-4910

Internet <http://www.protocol.com>

Technical Services email:
solutions@protocol.com

Marketing Department email:
marketing@protocol.com

Protocol Medical Systems Ltd.:

NW Europe

Derby Service Centre
St. Georges House, Vernon Gate
Derby DE1 1UQ
United Kingdom

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Improve NIBP Accuracy with Smartcuf™


Software Versions
2.0X-2.1X (NOTE: SMARTCUF button only available for 2.1X)

NIBP measurements can be adversely affected by many factor such as cardiac arrhythmias, sudden changes in blood pressure, body motions such as convulsions or shivering, bumping the cuff, vibration, vehicle motion, or weak pulses.

The patented Smartcuf software filtering technology greatly increases NIBP measurement accuracy in the presence of motion artifact or diminished pulses. Smartcuf synchronizes the NIBP reading with the R-wave of the patient's ECG to eliminate noise created by external stimuli such as patient motion or vibration. The monitor must perform ECG monitoring while using Smartcuf.

To enable the Smartcuf filter:


- Connect the ECG leads to the patient and perform ECG monitoring during NIBP.
- From the Main Menu, press **NIBP** to display the NIBP Menu and set Smartcuf to ON.

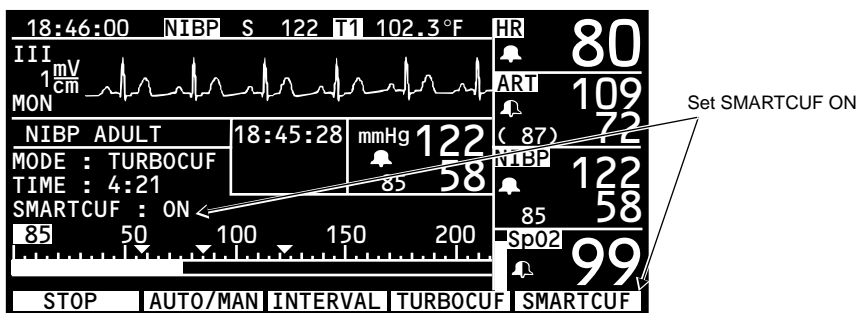English (cont.)

If artifact is so severe while Smartcuf is enabled that it affects the accuracy of an NIBP measurement, that measurement is marked with a special symbol on the display and on printouts

This symbol indicates the NIBP reading was taken in the presence of high artifact while monitoring ECG with Smartcuf enabled. Artifact can affect accuracy.

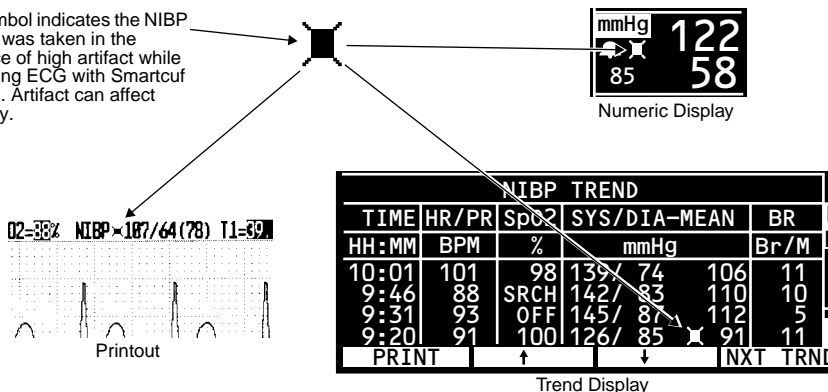There may be some situations where it is desirable to disable Smartcuf. This may include situations with very extreme motion artifact, certain types of arrhythmias, or other situations where it is not possible to obtain a good ECG signal. NIBP measurements can still be performed when Smartcuf is disabled.

To disable Smartcuf, from the Main Menu press **NIBP** to display the NIBP Menu and set Smartcuf to OFF.

NIBP Equipment Alert Message


Software Versions

2.1X

ARTIFACT PRESENT, MINIMIZE ARTIFACT (ERR# 15). The monitor has detected too much artifact to allow accurate readings. Take steps to reduce artifact. Position the patient's limb away from the body so the applied cuff is not in contact with the patient's body or any other object such as a bed rail. If the Smartcuf motion artifact filter is on, make sure that the ECG leads are properly connected to perform ECG monitoring during NIBP. If the Smartcuf motion artifact filter is off, consider turning it on (and connect ECG if not already connected).

NIBP Cuff Symbols

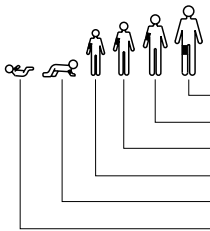Software Versions
2.0X-2.1X

 <p>NIBP cuff sizes:</p> <ul style="list-style-type: none"> Thigh Large adult Adult Small adult Child Infant 		<p>Apply the NIBP cuff as shown.</p>
		<p>Single-use only (not reusable).</p>

NIBP Display Default Settings

To select which vital sign waveforms are displayed, press **MAIN MENU**, **SETUP**, **WAVE SEL** to display the Wave Select window.

SETUP	WAVE SELECT	(85) 58
ECG : ON	RESP : OFF	PA 35/ 18
ART : ON	SpO2 : ON	MC02 BR
PA : OFF	NIBP : ON	35 12
C02 : OFF		mmHg
		SpO2 92
NEXT	ON/OFF	INSERV
		PREV MENU

To enable the display of a waveform or large NIBP numerics, select ON (ECG cannot be set to OFF). The monitor displays the first three active waveforms set to ON in the order of priority listed in the Wave Select window. If NIBP is ON and only one or two other waveforms are ON and active, the monitor displays large NIBP numerics in a waveform window.


The factory default settings for waveform display are:

Software version 2.03.00-2.1X	All waveforms and NIBP are ON.
Software version 2.00.00-2.02.00	All waveforms are ON. NIBP is OFF.

Press **NET OFF** to Disconnect from Acuity

Software Versions	Access to NET OFF
2.03.00-2.1X	NET OFF button on front panel NET OFF softkey (SETUP , ACUITY , NET OFF)
2.00.00-2.02.00	NET OFF softkey (SETUP , ACUITY , NET OFF)

When you want to disconnect the monitor from Acuity, be sure to use the **NET OFF** button as described below. The **NET OFF** button is a safety feature to help make sure the patient is not disconnected accidentally. If you simply disconnect the Acuity network cable from the monitor without using the **NET OFF** button, the monitor and Acuity both generate equipment alert messages that must be acknowledged by an operator.

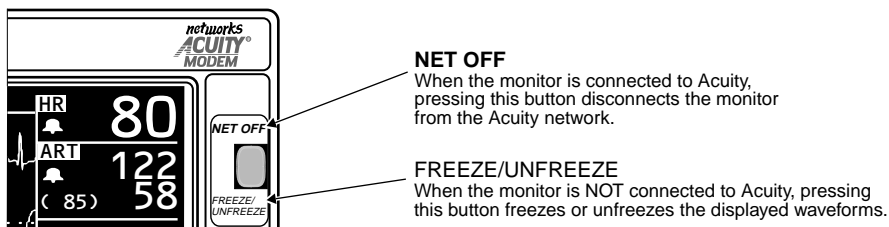1. To disconnect the Propaq Encore from the Acuity network, press the front panel **NET OFF** button.
2. Within 15 seconds, disconnect the Acuity network cable from either the Propaq Encore side panel or the bedside jack. If the patient will no longer be monitored with the Propaq Encore, turn off the monitor to erase trend information.

If you do not disconnect the Acuity network cable within about 30 seconds, the monitor and Acuity attempt to reconnect and prompt you to confirm the patient ID.

Whenever the monitor is connected to Acuity, the **NET OFF** function overrides the **FREEZE/UNFREEZE** function for the **NET OFF** front panel button. Therefore, you cannot freeze the displayed waveforms when the monitor is connected to Acuity. Also, you cannot initiate a Freeze Print at Acuity when connected to Acuity. When the monitor is not connected to Acuity, pressing the **FREEZE/UNFREEZE** button freezes or unfreezes the displayed waveforms.

Capnography (CO₂)

Software Versions
2.0X-2.1X

General CO₂ Specifications (Mainstream CO₂ and Sidestream CO₂)

Characteristic	Specification
CO₂ Display	
Accuracy	Mainstream¹: 0-30 mmHg, ± 3 mmHg 31-99 mmHg, $\pm 10\%$ of value Sidestream²: 0-30 mmHg, ± 3 mmHg 31-99 mmHg, $\pm 10\%$ of value
Breath Rate Display	
Accuracy	± 1 breaths/min or $\pm 5\%$, whichever is greater ³
CO₂ Performance	
Specification	Per ISO 9918:1993 (E) / EN 864:1996

1. Based on these airway conditions: sensor temperature = 42°C, airway adapter temperature = 33°C, water vapor pressure = 38 mmHg; standard gas mixture = CO₂ in balance air, fully hydrated at 33°C; barometric pressure = 760 mmHg and flow = 60 ml/min.
2. Based on the following additional airway conditions: Sample line = 7 ft, 0.055 in ID (2.13 m, 1.4 mm ID); Sample flow rate = 175 ml/min; Protocol watertrap (new/unused); Respiratory rate ≤ 50 bpm, stable to ± 3 breaths/min; Inspired/Expired time ratio = 1:2; Barometric pressure = 760 mmHg.
3. For Sidestream CO₂, this applies only for BR ≤ 50 .

Sidestream CO₂ Monitoring

Warning

If the Sidestream CO₂ option is connected to a ventilatory circuit, be sure to adjust appropriate ventilator or anesthesia system settings to compensate for the sampling flow volume (90 or 175 ml/min) that is aspirated from the ventilatory circuit by the Sidestream CO₂ option.

INCO₂ Numeric Display

The monitor only displays the numeric value for inspired CO₂ (INCO₂) if it is in alarm or if it is greater than or equal to 7.5 mmHg (or ≥ 1 kPa or 1%). Refer to the CO₂ specifications in the *Propaq Encore Reference Guide* for more information.

Time/Day Settings and Trends

Software Versions
2.0X-2.1X

Warning

Changing the hour/minute/second setting for the monitor in the Time/Day window can cause the monitor to erase previously stored patient trend data.

When you change the hour/minute/second setting for the monitor in the Time/Day window, the monitor deletes any patient trend data that is older than five hours for non-NIBP trends or older than eight hours for NIBP trends according to the new clock setting.

However, if the monitor has not yet stored the full capacity of trends and you change the hour/minute/second setting to a time that is within the stored trend period, previously stored trends are not erased.

Changing the day, month, or year setting does not affect the stored patient trends.

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